

## 2.2.7 Installing the AGSS Assembly (Optional)

1. Install two bracket dowel pins on the guide rail of the AGSS bracket, as shown in Figure A.
2. Use two M4×12 countersunk head screws to install the guide rail of the AGSS bracket on the left plate of the trolley column, as shown in Figure B.
3. Put the AGSS assembly on the AGSS bracket guide rail on the side plate of the trolley along the guide rail slot and tighten the knob, as shown in Figure C.
4. Install the AGSS transfer hoses at the exhaust gas outlet of the main unit and the AGSS assembly outlet, as shown in Figure D.

Figure A




Bracket dowel pin

Figure B



Guide rail of the AGSS bracket

Figure C

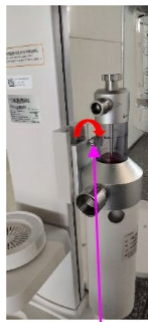

Knob

Figure D



AGSS transfer hose

May need to order a ¼" purple vacuum hose with that nut & nipple connector

